



Medical Carts

Comprehensive Product Catalog



The Avalo[®] Difference

Meeting the Needs for a Variety of Care Applications

With an emphasis on design, functionality, and organization, **Avalo Medical Carts** are built to support the way clinicians work. Avalo Carts can be deployed in hospitals, surgery centers, physicians' offices, assisted living facilities, and more. To support different care needs, Avalo carts can be custom configured to meet specific facility applications.

1. **Durable, high-impact panels**, built for long-term reliable performance. Stain, scratch and dent resistant.
2. **Medications and supplies are secured** with a selection of innovative locking systems including key, keyless or quick-access breakaway locking.
3. **Flexible divider systems and removable drawer trays** provide optimal visibility, grouping and organization of medications and supplies.
4. **Seamless drawers, smooth panels and work surface** promote simple cleaning for infection prevention.
5. **Accessory options** position equipment and supplies close at hand for easy access and improved efficiency, and keep the workspace clear.
6. **Three drawer depths** provide ample capacity to store and organize medications and supplies.
7. **A lightweight ergonomic design** and spread caster position provides stability control and ease of maneuverability.



Specifications

Standard Height



Intermediate Height



Compact Height



Drawer Specifications

	Height	Depth	Width
Three Inch	3" / 7 cm	17.5" / 45 cm	21.5" / 55 cm
Six Inch	6.25" / 16 cm	17.5" / 45 cm	21.5" / 55 cm
Ten Inch	10" / 25 cm	17.5" / 45 cm	21.5" / 55 cm

Cart Specifications

	Height	Depth	Width
Standard Height:	40.1" / 110 cm	24" / 61 cm	31" / 79 cm
Intermediate Height:	39.5" / 100 cm	24" / 61 cm	31" / 79 cm
Compact Height:	36" / 91 cm	24" / 61 cm	31" / 79 cm

Avalo® Medical Carts are Certified



RoHS Restriction of Hazardous Substances



ISO Quality Management System Standard



ISO Quality Management System Standard



UL Meets UL Safety Requirements



TUV Recognized for Quality, Safety and Sustainability

Avalo[®] SECURITY OPTIONS

LOCK SYSTEM CHOICES



Key Access with Manual Relocking (standard)

Featuring the Best[®] Core Removable Lock System, this basic yet proven lock system provides your medical cart advanced security and the flexibility to change keys simply and easily.



Breakaway Lock

Avalo's unique breakaway locking handle is standard on every model, offers the security you require and ensures quick access to code cart supplies



Keyless Access with Manual Relocking

Simply enter a 4-digit access code, push down the lock lever and the cart is ready for access. To secure the medical cart, simply push in the lock lever. Clear visual indication of lock status is apparent from a distance.



AutoLock[™] Automatic Relocking System

Automatically secures a motionless unlocked medical cart with a user-defined relock time delay from 1-99 minutes. AutoLock permits up to 5,000 individual user codes to further enhance cart access tracking and security. SLA battery operated.

AutoLock security upgrades available.

AUTOLOCK[™] UPGRADES



Proximity Access

An advanced proximity card reader can be selected as an additional security access requirement



Standard Host Software

Using Standard Host Software authorized users can add, change or delete an access code for specific users.

Build Your Custom Cart in 4 Steps

If specific drawer configurations or accessories are required to meet unique care applications, Avalo carts can be fully configured in a series of simple steps. Please contact your Capsa Representative for help configuring a cart to meet your needs.

STEP 1

Select your cart height & drawer configuration



Standard

Intermediate

Compact

CHOOSE DRAWER CONFIGURATION

Choose a drawer configuration combined with optional drawer dividers and trays to meet your storage needs.

STEP 2

Select your lock system (see page 3 for descriptions)



AutoLock™



Keyless Access with Manual Relocking



Key Access with Manual Relocking



Breakaway Lock

STEP 3

Select your color

Choose from our palette of colors:



Extreme Red




Extreme Blue



Extreme Green



Extreme Yellow

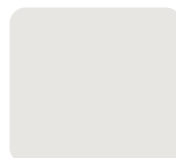

Light Crème



Dark Crème



Sample: Light Crème with Extreme Green Fillers

STEP 4

Select your accessories (see page 5)

Customize your cart with optional accessories to support the unique storage, organization and efficient workflow required.

Avalo[®] OPTIONS & ACCESSORIES



Cardiac Board



I.V. Pole



Equipment Shelf



Oxygen Tank Holder



Divided Drawer Tray



Drawer Dividers



Bulk Storage Bins



Tilt Bin Organizers



Side Organizer



Waste Container



Sharps/Glove Holder



Storage Container



Internal Lock Box



Writing Surface



Drawer Tag Holder



Tracking Caster



Accessory Bridge



Clear Top Mat



Plug Outlet Strip
UL1363A Rated



Computer Arm

M-Series TREATMENT CART

The M3-TC Treatment Cart offers a large storage capacity with easy maneuverability in a lightweight frame. Smart design provides plenty of room for supplies, and drawers are customizable for optimal organization. The M3-TC improves workflow and medical treatment preparation.

Features:

- Large capacity in a combination of shallow, medium, and deep drawers
- Integrated top handle and full-width drawer handles
- Reliable and easily changeable lock mechanism
- Durable design for long-term value
- 10 year limited warranty



Drawer Specifications

	Height	Depth	Width
Four Inch	4.25" / 11 cm	18.5" / 47 cm	22" / 56 cm
Six Inch	6.18" / 16 cm	18.5" / 47 cm	22" / 56 cm
Nine Inch	9.18" / 23 cm	18.5" / 47 cm	22" / 56 cm

Cart Specifications

	Height	Depth	Width
Cart Base	39.4"	26.75"	29.5"

Vintage Encore TREATMENT CART



The Vintage Encore Treatment Cart combines furniture styling with lightweight construction and impressive storage capacity. Smart design provides plenty of room for supplies, and drawers are customizable for optimal organization.

Features:

- Elegant styling to suit any facility's decor, in three furniture-style finishes
- Large capacity in a combination of shallow, medium, and deep soft-close drawers
- Simple-to-push casters
- No-bounce extended work surface
- Durable design for long-term value
- 10 year limited warranty

Available Finishes



Drawer Specifications

	Height	Depth	Width
Four Inch	4.25" / 10.8 cm	18.5" / 47 cm	23" / 58 cm
Six Inch	6" / 15.2 cm	18.5" / 47 cm	23" / 58 cm
Nine Inch	9" / 22.9 cm	18.5" / 47 cm	23" / 58 cm

Cart Specifications

	Height	Depth	Width
Cart Base	40.1"	25.2"	28"

Additional Information

Superior Products. Superior Service.

Contact Capsa Healthcare for more information about our wide selection of Medical Cart solutions that promote orderly storage and organization where it matters most. Adaptable to fit any application, Capsa's medical cart solutions are built to meet the unique needs of facilities worldwide.



CERTIFIED ISO
27001:2022



CERTIFIED ISO
9001:2015



CONTACT US

Contact us today for more information or to schedule a personal on-site demonstration.

800.437.6633

International: +1.614.864.9966

Email: Info@CapsaHealthcare.com

CapsaHealthcare.com